


Starting your personal therapy session


Note that the down arrow indicates more information on the next screen; press the DOWN arrow to scroll and view more details.

The AREA MENU will be presented at the start of each new area requiring treatment – only once if an individual area was selected. each treatment area can either be 1.STARTED, 2. SKIPPED or 3.RESTARTED. The system will automatically move to the next area in succession.

Now to commence treatment the display will ask you to move the device over the treatment area

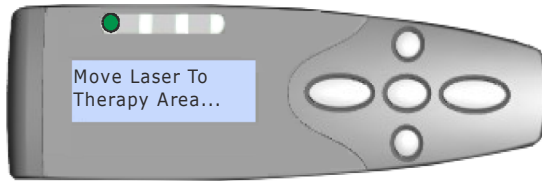

Once a surface is detected the lamp will turn on and the treatment will be administered. The RED LED will also be activated to indicate the Laser is active. The timer will display with the time REMAINING for this treatment session. A session can be STOPPED, PAUSED or SKIPPED by pressing ENTER during the session.

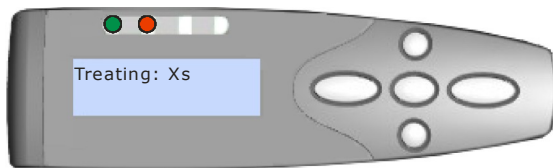

Congratulation! You should be now receiving treatment for your condition, and on the road to progress. Always remember to wear the safety goggles (as supplied) during all phototherapy procedures.

Important : Never repeat treatment more than once in a 24 hour period, unless under direct medical supervision.

Refer to instructions Manual or alternately contact your specialist if in doubt.

Safety

Your Dermaray Laser is a powerful and robust hand held Laser Device, As a Class IIIb laser, It must only be used in accordance with the Instructions as detailed within the User—Instruction Manual.” Incorrect and improper usage of this device may be dangerous and will void your warranty.

Caution When Using Laser Technology

Low Level Lasers produce no heat – any heat felt is generated by the actual electronic components. There is no evidence to suggest over exposure is possible or side effects exist. However, it is recommended to abide by the suggested exposure times. Longer exposure times may not produce better results from treatment. If symptoms persist, please see your Dermaray care professional.

Eye and accidental eye over exposure should be avoided with laser therapy. Laser light may cause permanent damage to your eyes. Damage may be cumulative, so repetitive accidental exposure may cause long term problems. So Please ensure that you use your Safety goggles for all treatment sessions.



For your convenience you can view some on-line instructions of how to use your Dermaray Laser by visiting the following web site links

Charging the battery

<http://www.youtube.com/watch?v=-Kmoz10sj2g>

Completing the Setup Wizard

www.youtube.com/watch?v=7w-xMKVVpdU

Instructions on How to Start and Use

www.youtube.com/?v=IEDHknYmxol

Total Control Pty. Ltd.
3/6 Barnett Court Morley
WA 6062 Australia
Web www.dermaray.com
All rights reserved.

Quick Start Instructions Guide



Dermaray Laser Comb 100mW Model DM-100



All Dermaray Branded products are designed and manufactured by Total Control Pty. Ltd. patents pending
All rights reserved 2010

Setup

Your Dermaray - Laser comb is different from other laser phototherapy devices. It is battery operated, contains an inbuilt computer including a real time clock and calendar, is fully programmable, very versatile and comes with a universal power supply for charging the batteries.

Your Dermaray - Laser comb intelligently manages internal power and operates the laser modules at variable frequency for maximum efficiency. This sustains the energy output while also increasing the life span of the laser modules resulting in more effective treatment for you.

Prior to commencing treatment or turning on the device, you should fully charge the battery in your Dermaray - Laser for a full 8 hours. (refer to instructions section 5) Download your copy from www.dermaray.com

How to use these instructions


The quick-start guide will enable you to become familiarized with your new Dermaray Laser Comb, Before attempting to use or charge the device – You must download your **free copy of the complete Instruction Manual from the Dermaray Web site**. Follow the link to the download section enter your details as requested, then download your copy.

Spend some time reading the manual, Follow the full instructions and then proceed with setup as required. If you are using a Laser Hair device for the first time. Please follow the suggested process as listed below;

- **Read** the user Instructions Manual
- **Charge** batteries as described in the instructions
- **Configure** - Laser for initial use
- **Refer** to the Instructions guide as required
- **Assistance**— contact your Dermaray consultant
- **Technical Support**— contact Dermaray

Your Dermaray Laser Comb may be pre-configured in different ways for different treatment purposes, including both strength and wavelength of the laser modules. Red and Blue laser are visible laser light, where as infrared lasers are invisible light, but are still harmful to the eye. If you Dermaray Laser is pre-configured to your treatment plan, do not adjust or tamper with it in anyway, always refer back to you Dermaray consultant.

Laser Comb features

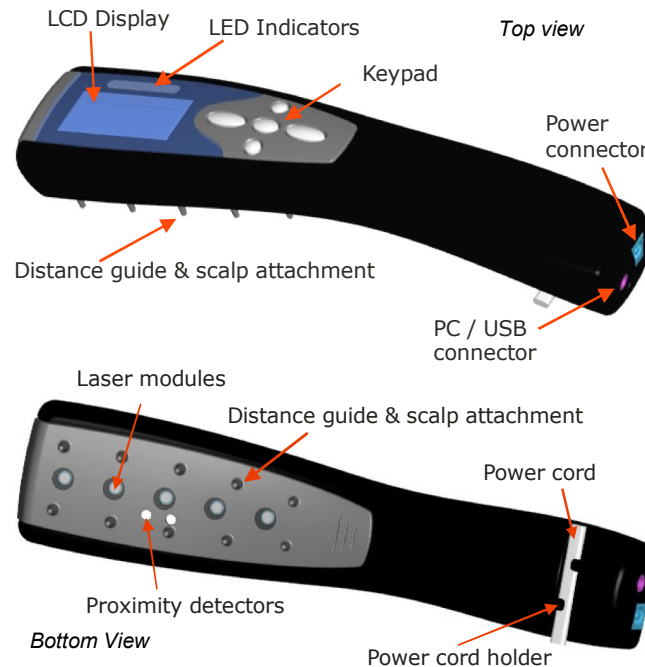

Application for optimal results

When using your Dermaray Laser it should be held against the treatment area and moved slowly around the surface for the allocated time (as per instructions from your therapist). For stubborn areas - the Laser can be held in one place for several minutes or for the duration of treatment if required. Moving - the Laser too quickly may result in a warning alert and deactivation of the lasers.


Laser Set-up Wizard

The WIZARD is an automated setup procedure; it will prompt and guide you during the setup phase, and request you to enter your patient treatment information. It is suggested to make some notes prior to starting.. The Setup Wizard must be completed within 25 minutes otherwise you will need to start again. (For detailed Information on how to enter the SET-UP patient information stage, refer to the User Instruction Manual and follow the flow chart guide)

If your Laser Comb is already pre-configured, then its time to start therapy


Starting your personal therapy session

When starting therapy setup you main screen look like the window below.

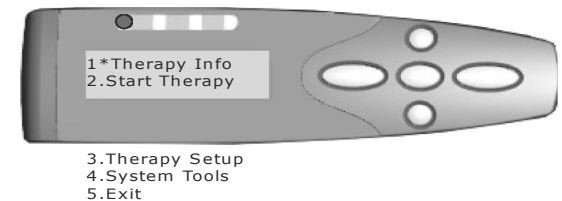

From the MAIN MENU, Select

1. Therapy-Info will display a list of patient names.
2. Select your name if indicated and press enter
3. Now select the area for treatment from your list and press enter.

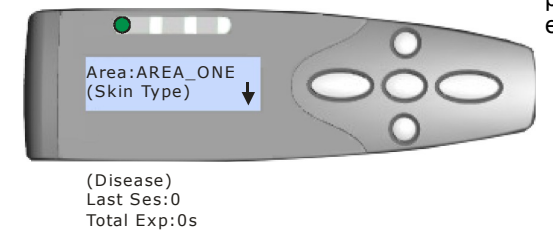

Information for each treatment area will be displayed for you. examples as per window.