



Starting your personal therapy session

At the TREATMENT AREA menu select 1.All Areas to treat all pre-programmed areas or select an individual area to perform treatment on.

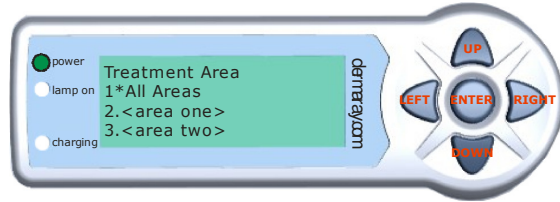

After making a selection at the TREATMENT AREA MENU above, the AREA MENU will be presented at the start of each new area requiring treatment – only once if an individual area was selected. Each Treatment Area can either be 1.STARTED, 2. SKIPPED or 3.RESTARTED. The system will automatically move to the next area in succession

At the completion of inputting all treatment areas move Lamp to the area for treatment.

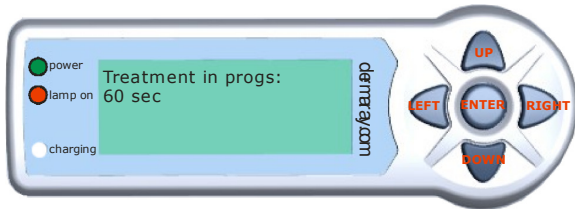

Once a surface is detected the lamp will turn on and the treatment will be administered. The RED LED will also be activated to indicate the UV lamp is active. The timer will display with the time REMAINING for this treatment session. A session can be STOPPED, PAUSED or SKIPPED by pressing ENTER during the session

Important : Never repeat treatment more than once in a 24 hour period, unless under direct medical supervision.

Refer to instructions Manual or alternately contact Dermaray or your supplier if in doubt.

Safety

Dermaray - UV is a powerful hand held Ultra Violet (UV) lamp, it is a Class IIa ionizing radiation emission device. It must only be used in accordance with the User's Instruction Manual." Incorrect usage of this device may be dangerous to your health, and void your Warranty Status .

Caution When Using UV Technology

Sunburn from over exposure should be avoided with all phototherapy treatments. Similarly avoid excessive direct natural sunlight before, after and between treatments. If after treatment your skin appears red or pink, feels excessively itchy, hot or prickly discontinue use, until your skin has fully recovered. If symptoms persist, please see consult a medical professional. If burning results shorten over the next several treatments then consider extra recovery time between treatments. Phototherapy programs are a suggestion and should be used as a guide only.

Eye and accidental eye over exposure should be avoided with all phototherapy treatments. So please ensure that you use your Safety goggles (as supplied) for all treatment sessions. (refer to user manual for detailed instructions).



For you convenience you can view on-line instructions of how to use your Dermaray UV lamp by visiting the following web site links

Instructions on How to Start and Use
www.youtube.com/?v=IEDHknYmxol

Charging the battery
<http://www.youtube.com/watch?v=-Kmoz10sj2g>

Completing the Setup Wizard
www.youtube.com/watch?v=7w-xMKVVpdU



Total Control Pty. Ltd.
3/6 Barnett Court Morley
WA 6062 Australia
Web www.dermaray.com
All rights reserved.

Quick Start Instructions Guide

UV Lamp Model NB-311



**Portable Ultraviolet Light Therapy
for treatment of
Psoriasis, Vitiligo, Dermatitis
and other skin conditions**

All Dermaray products are designed and manufactured by
Total Control Pty. Ltd. patents pending
All rights reserved 2010

Setup

Your **Dermaray® UV Lamp** is a premium quality and “**Best of Breed**”, **medical device** that will deliver an optimum amount of UV therapy, to maximize the effectiveness of your treatment as advised or prescribed by your medical Specialist.

Your Dermaray® UV Lamp has numerous Intelligent and superior built-in features and safeguards, that no other phototherapy device can provide; It includes; an on-board **computer**, which will create your individual treatment plan, then, Monitor, Regulate, Control, Prompt, Record and Log, your usage. The **proximity detector**/sensors assist the user by alerting you of incorrect operation, and will automatically turn the device off, if it senses no activity, eliminating over- exposure, and possible burning. The ultra long life **re-chargeable** batteries ensures continuous and safe treatment sessions, without being directly connected to mains power, which can impede the users safety, as well as pose a danger; and the **Universal Power Supply** means that you can use your device anywhere in the world, (you will need only to acquire the correct power cord for your new country of destination)

Prior to commencing treatment or turning on the device, you should fully charge the battery in your Dermaray - UV for a full 6 hours, (refer to instructions manual). Download your copy now from **www.dermaray.com**.

How to use these instructions

The quick-start guide will enable you to familiarize yourself with the Dermaray UV Lamp, Spend some time reading the User Instruction manual, which you must download from www.dermaray.com now. Follow the full instructions and then proceed with setup as required. Please follow the suggested process as listed below;

- ⇒ **Read** the user Instructions Manual
- Charge** - batteries as described in the instructions
- Configure** - UV Lamp for initial use
- Refer** - to the Instructions guide as required
- Assistance** - contact Dermaray or your distributor
- Technical Support** - contact Dermaray or distributor

Operations

Your Dermaray lamp is a UV lamp appliance for the treatment of skin diseases such as psoriasis, vitiglio dermatitis etc. and hair loss. It must be used in strict accordance with your user instruction Manual. The UV Lamp may be configured in different ways for different treatment purposes, If your Dermaray UV is pre-configured by your medical specialist, do not adjust or tamper with it in anyway, always refer back to your place of purchase.

Features of UV Lamp



What's in the box

Your Dermaray UV Pack is supplied with the following:-

- ⇒ A Dermaray UV lamp
- Universal Mains power supply / battery charger
- UV lamp safety goggles compact (patient)
- UV lamp safety goggles large (patient / helper)
- Power cord
- Clear plastic distance guide and scalp attachment
- Quick-Start user guide


UV Lamp Set-up Wizard

The WIZARD is an automated setup procedure; it will prompt and guide you during the setup phase, and request you to enter your patient treatment information. It is suggested to make some notes prior to starting.. The Setup Wizard must be completed within 25 minutes otherwise you will need to start again. (For detailed Information on how to enter the SET-UP patient information stage, refer to the User Instruction Manual and follow the flow chart guide)

If your UV Lamp is already pre-configured, then its time to start therapy

Starting your personal therapy session

When starting treatment setup, you main screen look like the window below. WAKE Dermaray UV lamp, by pressing and holding the ENTER key for 3 seconds. Then at the WELCOME screen, press and hold the RIGHT key for 3 second duration.

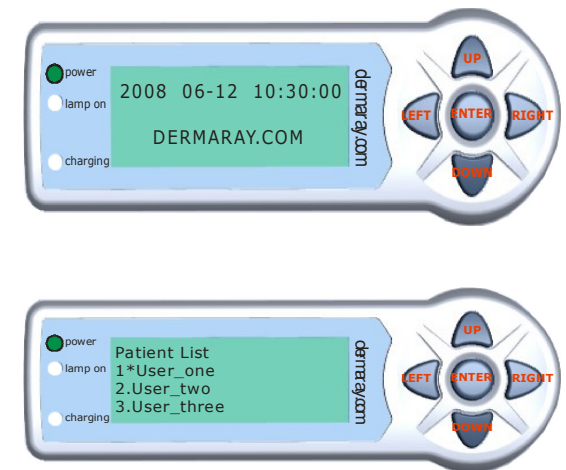

From the MENU, Select

1. Patient Name (select your name and press enter)
2. Now select the area for treatment from your list and press enter.